

# Pathfinder Plus™

Surgeon Controlled  
Urologic and  
Hysteroscopic  
Irrigation System



# Surgeon Control.

The surgical team benefits when the surgeon has control of irrigation. The surgeon knows exactly when and how much bolus irrigation is required. The surgical assistant can then focus on other critical tasks.

Pathfinder Plus is the elegantly simple irrigation device that provides this control. The result is greater accuracy and enhanced visualization during delicate endoscopic procedures.

## Improved Visualization


Ensure a clear view by quickly clearing debris from the surgical field. Laser beam scatter is also reduced, providing accurate delivery of energy.

## Selective Irrigation

Pathfinder Plus provides a slow continuous flush which is controlled by saline bag height and/or pressure cuff. In addition, bulb compression provides a precise flush bolus — a gentle squeeze is effective for clearing large amounts of debris from the field, while 'fluttered' pulses can be effective at altering stone position for improved surgical management.

## Universally Compatible

Standard Luer-lock fittings attach directly to any cystoscope, hysteroscope, or any other endoscope where the surgeon requires or benefits from improved results during delicate procedures.

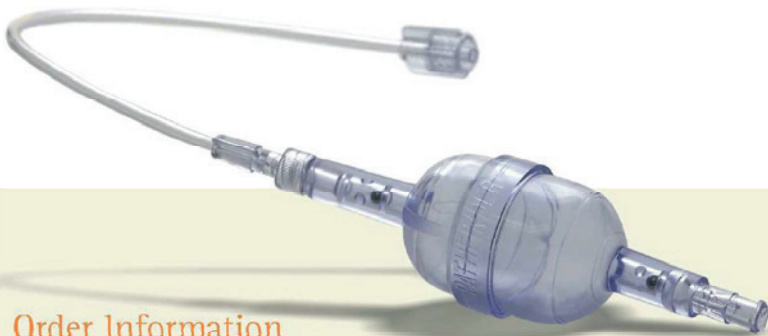

## Order Information

Pathfinder Plus may be ordered by:

- ◆ calling Customer Service at 800-533-4984
- ◆ faxing your order to 801-566-2062
- ◆ ordering online at [www.order.utahmed.com](http://www.order.utahmed.com)

Order No	Description	Package
PA701	Pathfinder Plus Endoscopic Irrigation System	10 / box


# COMMUNICATION

*Pathfinder Plus eliminates communication problems between surgeon and assistant by providing the surgeon a precise and unobtrusive means of irrigation control.*



# CONTROL

*Pathfinder Plus bulb fits in scope hand, preserving dominant hand for instrumentation.*

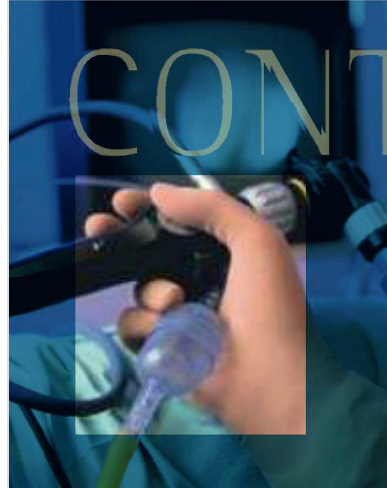

# CONVENIENCE

*Setup of Pathfinder Plus is simple, eliminating the time and frustration of setting up a complex system of tubing, footpedals, or syringes:*

- ◆ Attach Pathfinder Plus tubing to endoscope
- ◆ Attach Pathfinder Plus bulb to irrigation
- ◆ Squeeze bulb to flush system.



**Surgical Device GmbH**  
Alte Steinhäuserstrasse 19  
CH-6330 Cham  
Telefon +41 (0)41 784 38 00  
Telefax +41 (0)41 784 38 09  
[info@surgicaldevice.ch](mailto:info@surgicaldevice.ch)  
[www.surgicaldevice.ch](http://www.surgicaldevice.ch)